

EXAMINATION CHAIR Model 2000-CH

RIGHTmed
S4OPTIK



S4OPTIK

FEATURES:

◆ effortless CRADLE TILT RECLINE

precision engineered Cradle Track System provides effortless recline with perfect counterbalance.

◆ convenient TILT RELEASE / CONTROL

one-handed release/control mechanism for precise recline control, accessible on both sides of the chair, for easy, convenient patient positioning.

◆ strong, yet silent POWER LIFT SYSTEM

powered by Swiss designed linear actuator for smooth and virtually silent lifting with a 500lb (227 kg) capacity. The chair seat raises to 31" (78.7cm) and lowers to 23" (58.4cm).

◆ intelligent ROTATION/LOCK

foot actuated Rotation Lock provides a firm positive lock. Accessible from both sides of the chair, our unique design eliminates operator bending or reaching.

◆ sensible HEIGHT CONTROL LOCATIONS

height adjustment controls accessible from the slit lamp arm, stand console (when paired with a S4OPTIK Instrument Stand), or movable foot switch.

◆ convenient FOOT CONTROL

robust foot switch provides raise/lower control for chair height adjustment.

◆ advanced HEADREST

single handed control for easy positioning of the Headrest, insuring head and neck support. Engineered for a positive, non-slip lock in any location.



SPECIFICATIONS

Measurements

Base	22" (55.8 cm) wide x 23.25" (59 cm) depth
Seat	20" (50.8 cm) width, lower limit 23", upper limit 31"
Footrest	11.75" (29.8 cm) wide x 8" (20.3 cm) depth.
Shipping Weight	276 lbs (125.19 kg)

TECHNOLOGY
VALUE
RELIABILITY




◆ functional ARMRESTS

upholstered for comfort, plus individually retractable for patient flexibility. Solid construction provides support for ingress/egress.




◆ robust FOOT REST

a large, broad and stable foot rest provides safe and easy ingress/egress. Easily able to withstand the full weight of the patient and conveniently retracts when not required.

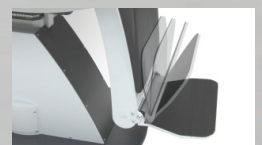

CHAIR UPRIGHT

CHAIR RECLINED

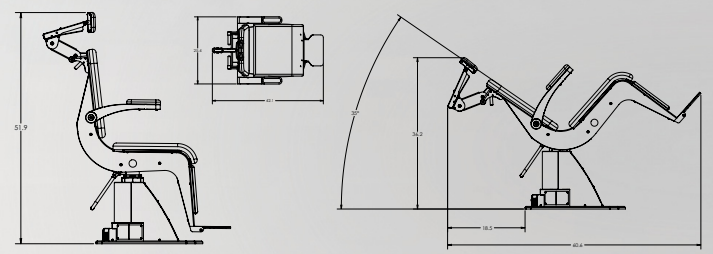

Distributed by:

with a policy of continuous improvement, we reserve the right to make changes to specifications without advance notice.

RIGHTmed
S4OPTIK

5325 Cleveland Street, Suite 303 Virginia Beach, VA 23462
1.888.224.6012 | www.s4optik.com